


**Dumont Medical Micro-Blunted
Atraumatic Tipped Forceps**

The sides of these forceps are distinctively ribbed to provide a better holding surface for maximum control of the instrument. The shanks are arched toward each other to transmit more holding strength to the points. The smooth inside surfaces at the tips will close flat to 5 mm with gentle pressure, increasing the total holding area considerably. The microscopically matched tips are gently blunted to prevent tissue puncture and are recommended as micro needle holders.

FORCEPS/TWEEZERS

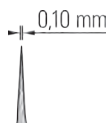

Medical #5

Standard
Straight

11 cm
Inox
0.1 x 0.06 mm

No. 11253-20



Medical #5A

Standard
Off-Center

11.5 cm
Inox
0.13 x 0.08 mm

No. 11253-21



Medical #7S

Standard
Curved

12 cm
Inox
0.17 x 0.1 mm

No. 11273-22




Medical #7

Standard
Curved

11.5 cm
Inox
0.17 x 0.1 mm

No. 11273-20

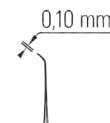

Medical #5/45

Standard
Angled 45°

11 cm
Inox
0.1 x 0.06 mm

No. 11253-25



Medical #5/90

Standard
Angled 90°

11 cm
Inox
0.1 x 0.06 mm

No. 11253-29