


**Dumont Medical Micro-Blunted  
Atraumatic Tipped Forceps**

The sides of these forceps are distinctively ribbed to provide a better holding surface for maximum control of the instrument. The shanks are arched toward each other to transmit more holding strength to the points. The smooth inside surfaces at the tips will close flat to 5 mm with gentle pressure, increasing the total holding area considerably. The microscopically matched tips are gently blunted to prevent tissue puncture and are recommended as micro needle holders.



FORCEPS/TWEEZERS

Medical #5

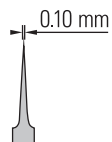

**Standard  
Straight**

11 cm  
0.1 x 0.06 mm  
Inox

No. 11253-20

Medical #5A




**Standard  
Off-Center**

11.5 cm  
0.13 x 0.08 mm  
Inox

No. 11253-21

Medical #7S

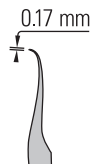

**Standard  
Curved**

12 cm  
0.17 x 0.1 mm  
Inox

No. 11273-22

Medical #7




**Standard  
Curved**

11.5 cm  
0.17 x 0.1 mm  
Inox

No. 11273-20

Medical #5/45

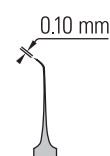

**Standard  
Angled 45°**

11 cm  
0.1 x 0.06 mm  
Inox

No. 11253-25

Medical #5/90



**Standard  
Angled 90°**

11 cm  
0.1 x 0.06 mm  
Inox

No. 11253-29