

Fine Forceps Repair Kit

Realign and repair fine forceps with this convenient kit. Bent forcep tips can be restraightened with the chain nose pliers. The sharpening stone has two grooves, 1.6 mm and 0.8 mm, to assist in resharpener the forcep tips. O-rings clamp the forceps shut to ensure that both shanks are equally aligned and sharpened. Instrument oil assists in producing a smooth sharp edge. Detailed instructions are included.



Kit includes:

- 1 Sharpening stone for fine forceps (No. 29008-22)
- Ten O-rings (11200-00)
- 1 Instrument oil, 50 ml bottle (No. 29055-00)
- 1 Chain nose pliers

No. 29000-00



O-Rings for Forceps

These rubber O-rings fit our forceps and keep them closed with variable degrees of pressure on the tips.

Pkg. of 10

No. 11200-00




Instrument Oil

Physiologically inert oil for lubricating surgical instruments, box locks, ratchets and any overlapping joints.

50 ml bottle

No. 29055-00

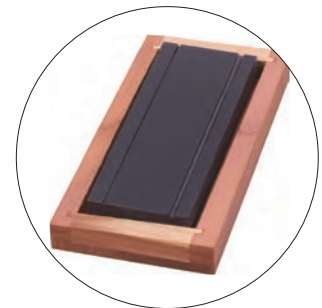

Sharpening Stone for Fine Forceps

This natural Arkansas stone assists in resharpener the tips of fine forceps. The stone has two grooves, 0.8 mm and 1.6 mm and comes in a cedar box with ten O-rings.

12.5 x 6.5 x 3 cm

No. 29008-22

INSTRUMENT CARE & ACCESSORIES